


Cue-Mate Procedure:

• Cue-Mate System and Features



Procedure Bulletin P1A

Wishbone

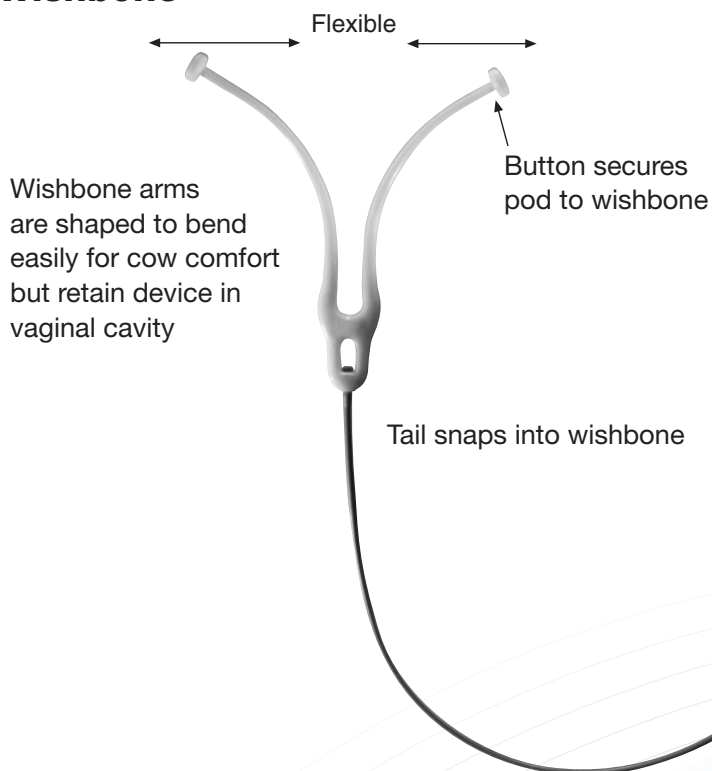

Progesterone Pod


Consists of biomedical silicon impregnated with progesterone

Cone provides soft round surface for comfort and safety


Gills are soft for comfort

Mucus flows over gills to dissolve and release progesterone

Pod design allows mucus to flow through ports, down inner surface maximising surface area for progesterone release

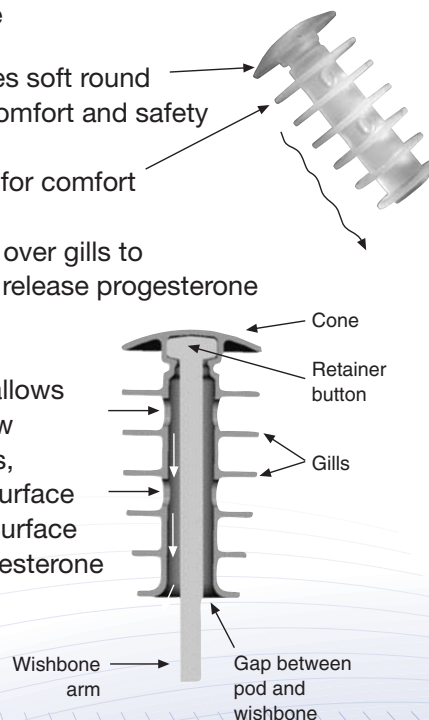

Cue-Mate components ready for assembly

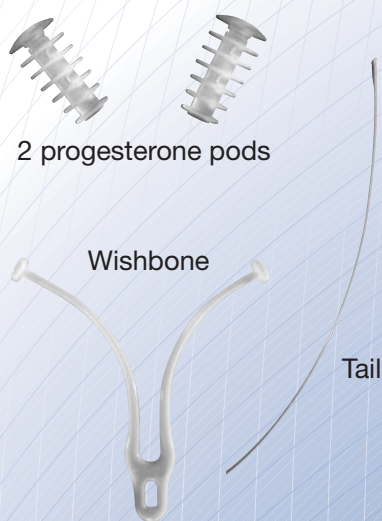

Cue-Mate assembled and ready to use



Cue-Mate Applicator (closed jaws)

